


VESTAKEEP® i5 G

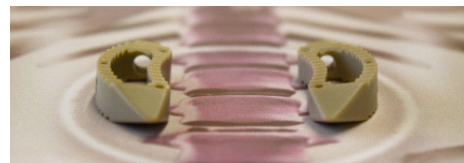

Implantable grade polyether ether ketone resin for permanent implants

VESTAKEEP® i5 G is a natural colored, very high viscosity polyether ether ketone (PEEK) that is especially designed for long term implantable medical devices.

Proven biocompatibility

The extra high purity and extended quality measures make VESTAKEEP® i-grade materials an excellent choice for permanent implants.

For VESTAKEEP® i5 G, biocompatibility has been tested following ISO 10993-1 recommendations for permanent tissue/bone contact and USP Class VI.

VESTAKEEP® i5 G complies to ASTM F2026 "Standard Specification for Polyetheretherketone (PEEK) Polymers for Surgical Implant Applications".

A summary of biocompatibility test results is available upon request.

VESTAKEEP® i5 G is supplied as cylindrical pellets.

Biocompatibility tests available for i5 G

Standard	Description
ISO 10993-3	Genotoxicity: Ames Test
ISO 10993-3	Genotoxicity: Chromosome aberration test
ISO 10993-3	Genotoxicity: Mouse Lymphoma test
ISO 10993-5	Cytotoxicity
ISO 10993-6	Test for local effects after Implantation in bone (90 days)
ISO 10993-10	Sensitization: Maximization test according to Magnusson and Kligman
ISO 10993-10	Irritation: Intracutaneous Reactivity
ISO 10993-11	Subchronic Systemic Toxicity
ISO 10993-12	GC/MS Fingerprint
USP Class VI	Acute Systemic Toxicity Intracutaneous Reactivity Muscle Implantation

Properties of VESTAKEEP® i5 G

		Test method	Unit	Value
Density	23°C / 50% r.h.	ISO 1183	g/cm ³	1.30
Water absorption	saturation	ISO 62	%	0.4
Moisture absorption	23°C / 50% r.h.	ISO 62	%	0.12
Tensile test	23°C / 50% r.h.	ISO 527-1 ISO 527-2		
Stress at yield			MPa	95
Strain at yield			%	5.0
Stress at break			MPa	> 70
Strain at break			%	> 40
Tensile modulus		ISO 527-1 ISO 527-2	MPa	3400
CHARPY impact strength	23°C / 50%r.h. -30°C	ISO 179/1eU	kJ/m ² kJ/m ²	N ¹⁾ N ¹⁾
CHARPY notched impact strength	23°C / 50% r.h. -30°C	ISO 179/1eA	kJ/m ² kJ/m ²	9 C ¹⁾ 8 C ¹⁾
Temperature of deflection under load				
Method A	1.8MPa	ISO 75-1	°C	150
Method B	0.45MPa	ISO 75-2	°C	205
Vicat softening temperature		ISO 306		
Method A	10N		°C	335
Method B	50N		°C	305
Differential Scanning Calorimetry (DSC)		ISO 11357		
Recrystallization temperature			°C	285
Glass transition temperature, 2 nd heating, onset			°C	145
Glass transition temperature, 2 nd heating, midpoint			°C	155
Melting temperature, 2 nd heating			°C	340
Melt volume-flow rate (MVR)	380°C/5kg	ISO 1133	cm ³ /10 min	7

¹⁾ C = Complete break, incl. hinge break H
N = No break

VESTAKEEP® is a registered trademark of Evonik Degussa GmbH

This information and any recommendations, technical or otherwise, are presented in good faith and believed to be correct as of the date prepared. Recipients of this information and recommendations must make their own determination as to its suitability for their purposes. In no event shall Evonik assume liability for damages or losses of any kind or nature that result from the use of or reliance upon this information and recommendations. EVONIK EXPRESSLY DISCLAIMS ANY REPRESENTATIONS AND WARRANTIES OF ANY KIND, WHETHER EXPRESS OR IMPLIED, AS TO THE ACCURACY, COMPLETENESS, NON-INFRINGEMENT, MERCHANTABILITY AND/OR FITNESS FOR A PARTICULAR PURPOSE (EVEN IF EVONIK IS AWARE OF SUCH PURPOSE) WITH RESPECT TO ANY INFORMATION AND RECOMMENDATIONS PROVIDED. Reference to any trade names used by other companies is neither a recommendation nor an endorsement of the corresponding product, and does not imply that similar products could not be used. Evonik reserves the right to make any changes to the information and/or recommendations at any time, without prior or subsequent notice.

Evonik Resource Efficiency GmbH

Kirschenallee, 64293 Darmstadt, Germany

Phone +49 2365 49-9227 e-mail evonik-hp@evonik.com